

EHBT-75

7-Diff Auto Hematology Analyzer

Deeper Clinical Insights from a Single Drop

- EHBT-75 is a 7-diff cell morphology analyzer designed to deliver intelligent blood diagnostics for human use
- It expands capabilities to include abnormal detection such as **NST, NSG, NSH, NLR, PLR, ALY, PAg and RET**, offering deeper clinical insights
- The system utilizes a single-use cartridge that can be stored at room temperature and is maintenance-free

Product Features

**AI-Powered Precision**

- Deep learning algorithms enable accurate cell identification and classification


**High-Resolution Morphology**

- Image-based detection with enhanced resolution captures clear cellular structures


**Liquid-Based Staining**

- Liquid-based cytology delivers richer color dimension


Results Visualization

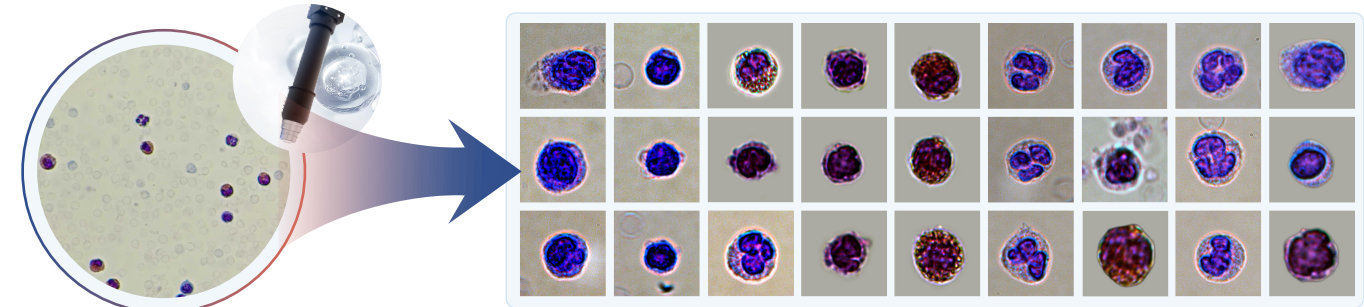

Fully Automated Sample Processing

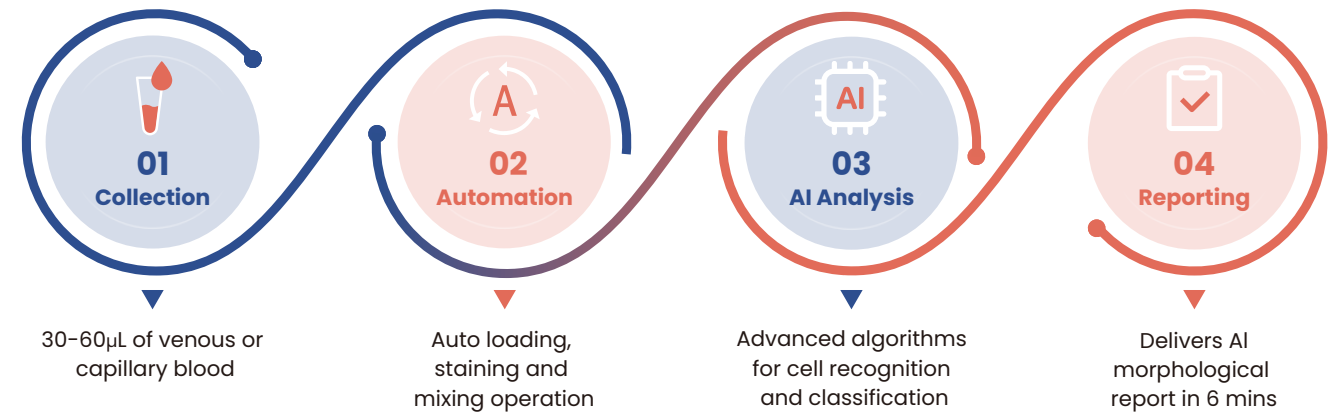

Single-Use Test Kit

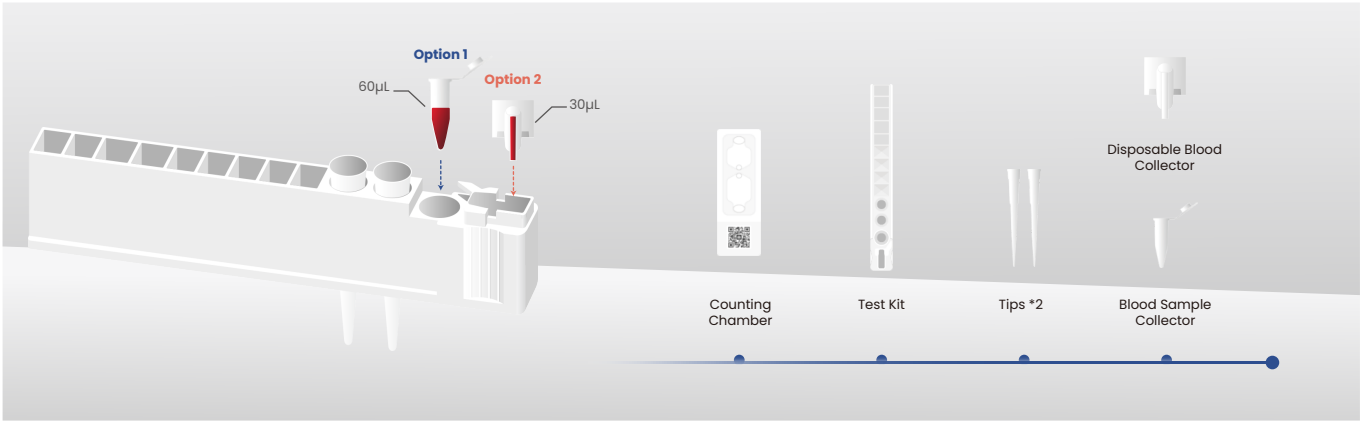

Specifications

Principle : Cell morphology for CBC, Photoelectric colorimetry for HGB

Parameters : WBC, NEU#, LYM#, MON#, EOS#, BAS#, NEU%, LYM%, MON%, EOS%, BAS%, RBC, HGB, HCT, MCV, MCH, MCHC, RDW-SD, RDW-CV, PLT, MPV, PDW, PCT, P-LCR, P-LCC, NST#, NST%, NSG#, NSG%, NSH#, NSH%, NLR, PLR, ALY#, PAg#, RET#, RET%

Sample Type : Capillary/venous whole blood

Sample Volume : 30-60µL for loading, 30µL for testing

Throughput : 10 samples/hour

Dimensions : 415mm × 203mm × 483mm (L × W × H)

Display : 7-inch screen

Printer : External thermal or A4 printer

Communication : 2 USB ports, Ethernet, WIFI, Bluetooth

Power Supply : AC 100-240V, 50/60Hz, 48VA

Ozelle

BEESKOWER STRABE 4 IN 15234 FRANKFURT, GERMANY

info@ozellepoct.com www.ozellepoct.com

CE

RPS/1653/2025

